VYNTUS™ QUICK GUIDE

JAEGER


Box Calibration Workflow Quick Guide

Vyntus™ BODY

1. Launch the Box Calibration program - two ways to access the Gas Calibration program



OR

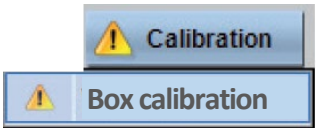

From the "Home Page"

From a Measurement or Calibration program


2. Click <**F7**> to activate the door magnets
   Close the body box door by pressing the door gently to the magnets




3. Click <**F1**> to "**Start**" the complete body box calibration, which will start automatically after 2 minutes stabilization time

   It includes the TAU verification and 3-frequency Box calibration



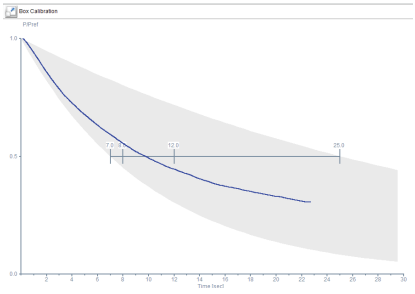

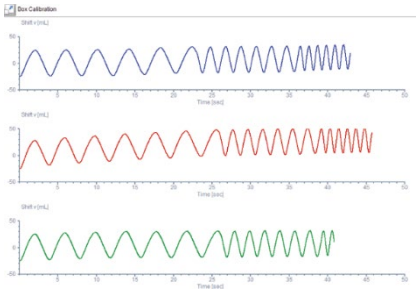

4. After the "**Box Calibration**" is completed, check if the results are within the limits = "OK"

   TAU between 8 and 12 seconds (acceptable up to 25)
   Gain between 0.75 – 1.5



5. Close the "**Box Calibration**" program - two ways to close the "**Box Calibration**" program



OR



Return to "Home Page"

Return to previous measurement program

jaegerrdx.com

For US Distribution © 2025 Jaeger Medical America, Inc. Jaeger and the Jaeger logo are trademarks or registered trademarks of Jaeger Medical America, Inc., or one of its affiliates. Please read the complete Instructions For Use (IFU) that come with the devices or follow the instructions on the product labelling. JAE-US-2500032.